Best Selling 808nm Diode Laser Hair Removal

Functions:

- 1. Hair Removal
- 2.Skin Rejuvenation
- 3.Skin care

Product exhibition:





Applications:

For fast, safe, painless and permanent hair removal on all 6 skin types, including dark skin. Suitable for any unwanted hairs on areas like face, arms, armpits, chest, back, bikini, legs...













Product descriptions:

Item	Best Selling 808nm <u>Diode Laser Hair</u> <u>Removal</u>
Type	Laser
Laser type	Semiconductor/freezing painless laser
Q-Switch	No
Style	Stationary
Laser instrument working mode	Pulse mode
Cooling Fin	240*2

Wave length	808nm+10nm
Light spot size	12mm*20mm
Pulse width	10-300ms
Maximum energy of a pulse	3J-140J
Pulses energy	140J+10J
Laser instrument power	600W
Input power	1200W
Working voltage	110V+ 10%,50Hz
OSD Language	Chinese, English
Insurance specification	Max 10 A
Size (Machine)	1180*500*385mm
Weight	≈51kg

Product detailed images:





Advantages:

- 1. Freezing point hair removal of semiconductor laser brings in **Germany micro channel cooling technology**. And **Germany imported bars**, confirm the shoting of the machine over 5,000,000 times and have enough energy.
- 2. **The professional filter system**, ensure the purity of the circulating water.
- 3. A precise constant temperature system, absolutely realize freezing point hair removal, fully guaranteed the sapphire skin protection system keep at the temperature of 4°C. And **strict temperature control**, if temperature over 35 degree, the machie will stop work, this can well protect the therapy handle safety, especiall for the bars and crystal filter.
- 4. 808nm semiconductor laser hand tool can realize quick and comfortable hair removal.
- 5. **12mm-20mm suitable big spot size**, it confirm the hair removal for any part of the body is flexible, such as lip hair removal, vagine hair removal, etc.
- 6. **Taiwan imported power supply**, this is the best power supply in the market, confirm the machine stable and energy enough.

Packaging&delivery:

Packaging	1.Styrofoam 2.Fumigation wooden case 3.Air box
Delivery	Based on your oredr quantity






Operation Principle:

- 1.The appropriate laser pulse time can ensure the accumulated energy to break the targeted tissue without damaging the surrounding tissues.
- 2. The appropriate energy density can provide enough energy output in a certain time to break the targeted tissues without damaging the surrounding tissues.
- 3. The appropriate epidermal protection measure can ensure enough damage of the targeted tissues without damaging the surrounding tissues.
- 4. The multi-pulse laser of this semiconductor laser hair

removal series chooses the mode of a relative low energy density to heat mascara to 75°C, and maintain this state with the hand tools sliding under 10Hz, then both the mascara and stem cell growth will become inactive, reaching the goal of permanent hair removal

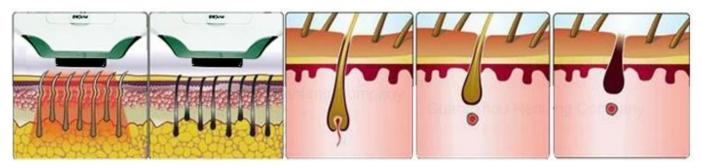

Contact: Racher Cheng

Tel&Whatsapp: 0086-15527106178

Skype: sales2@lanyuansupply.com

Email: sales2@lanyuansupply.com